

HEARING AID CLEANING GUIDE



DAILY




It is a good idea to wipe down the body of the hearing aid with the lint-free cloth provided in your kit. This is to remove any sweat or oils from your skin that have accumulated throughout the day. It is also a good idea to check the earpiece to see if there is any wax that needs wiped off. This is most commonly done before placing the hearing aids in the charger.

WEEKLY



We recommend using the brush provided in your kit to help loosen and remove dust or other small particles that could block the microphone ports. **DO NOT INSERT THE BRUSH INTO ANY OPENINGS ON THE HEARING AID.** Brush across the openings circled above plus an additional port on the opposite side of the hearing aid.

MONTHLY

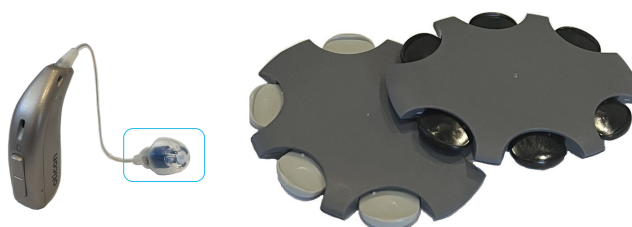

Once a month, check the wax trap to see if it needs replaced. Please see the wax trap handout for more information on where to locate the wax trap and how to change it.

WHEN TO MAKE AN APPOINTMENT

To keep your hearing aid working well and lasting longer, we recommend routine cleanings by your audiologist every 4-6 months. If you have tried all these cleaning steps and your hearing aids is still not working, call our office at (480) 331-6521 to make an appointment. Appointments are subject to office visit fees if you are outside your service plan period.